

This Page Is Inserted by IFW Operations
and is not a part of the Official Record

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images may include (but are not limited to):

- BLACK BORDERS
- TEXT CUT OFF AT TOP, BOTTOM OR SIDES
- FADED TEXT
- ILLEGIBLE TEXT
- SKEWED/SLANTED IMAGES
- COLORED PHOTOS
- BLACK OR VERY BLACK AND WHITE DARK PHOTOS
- GRAY SCALE DOCUMENTS

IMAGES ARE BEST AVAILABLE COPY.

**As rescanning documents *will not* correct images,
please do not report the images to the
Image Problem Mailbox.**

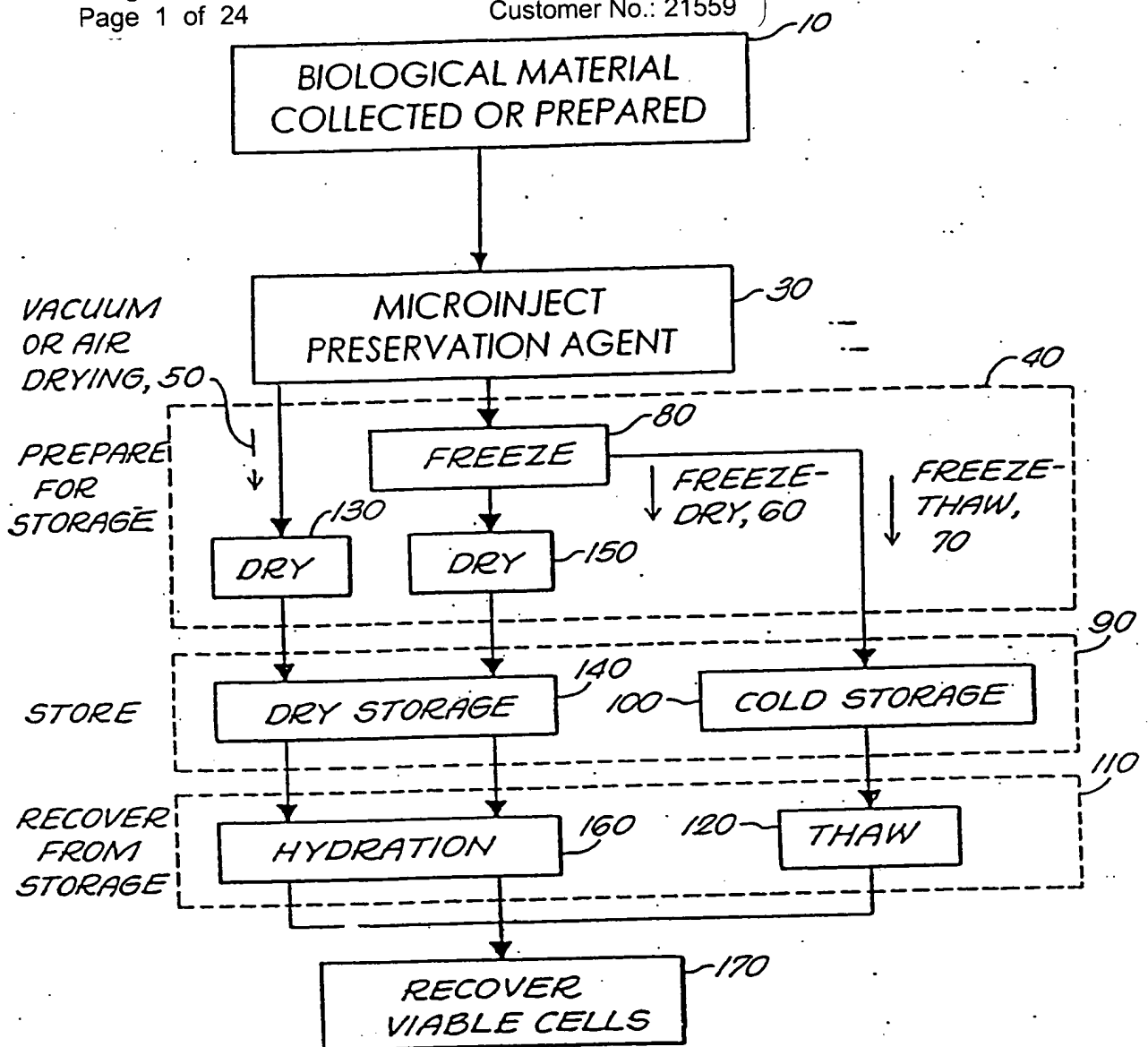


Figure 1



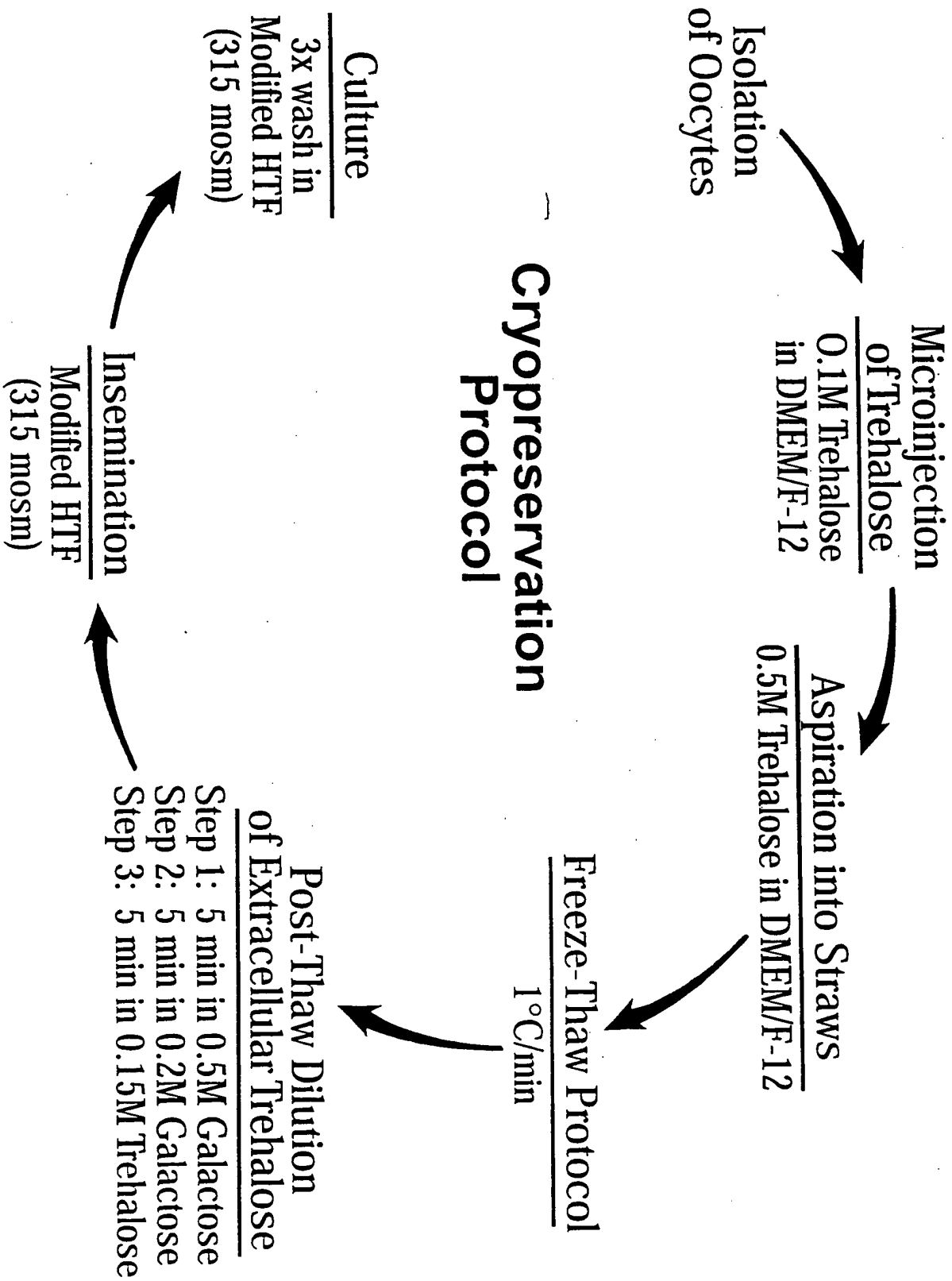


Figure 2

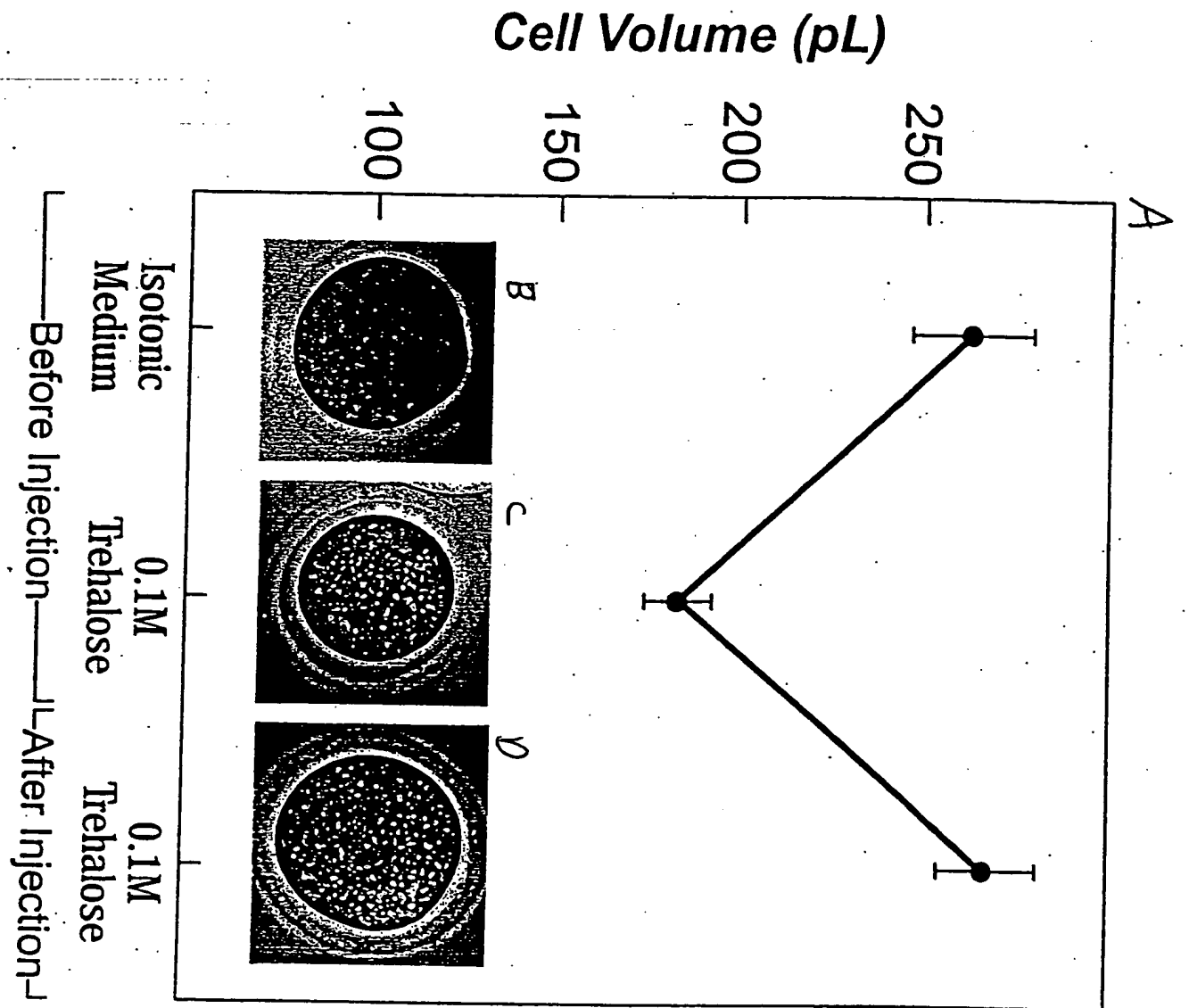

Figures 3A – 3D

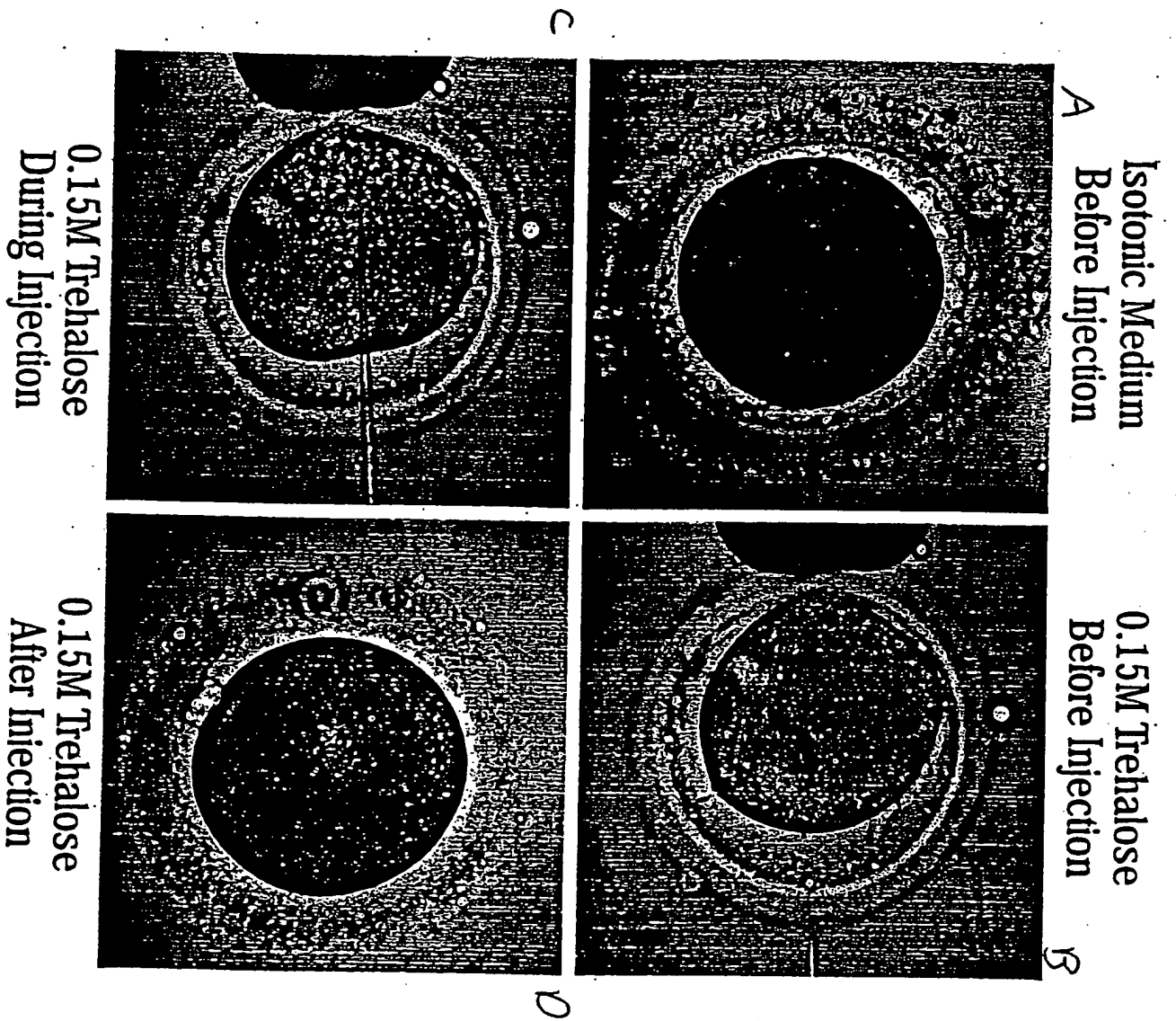

Figures 4A - 4D



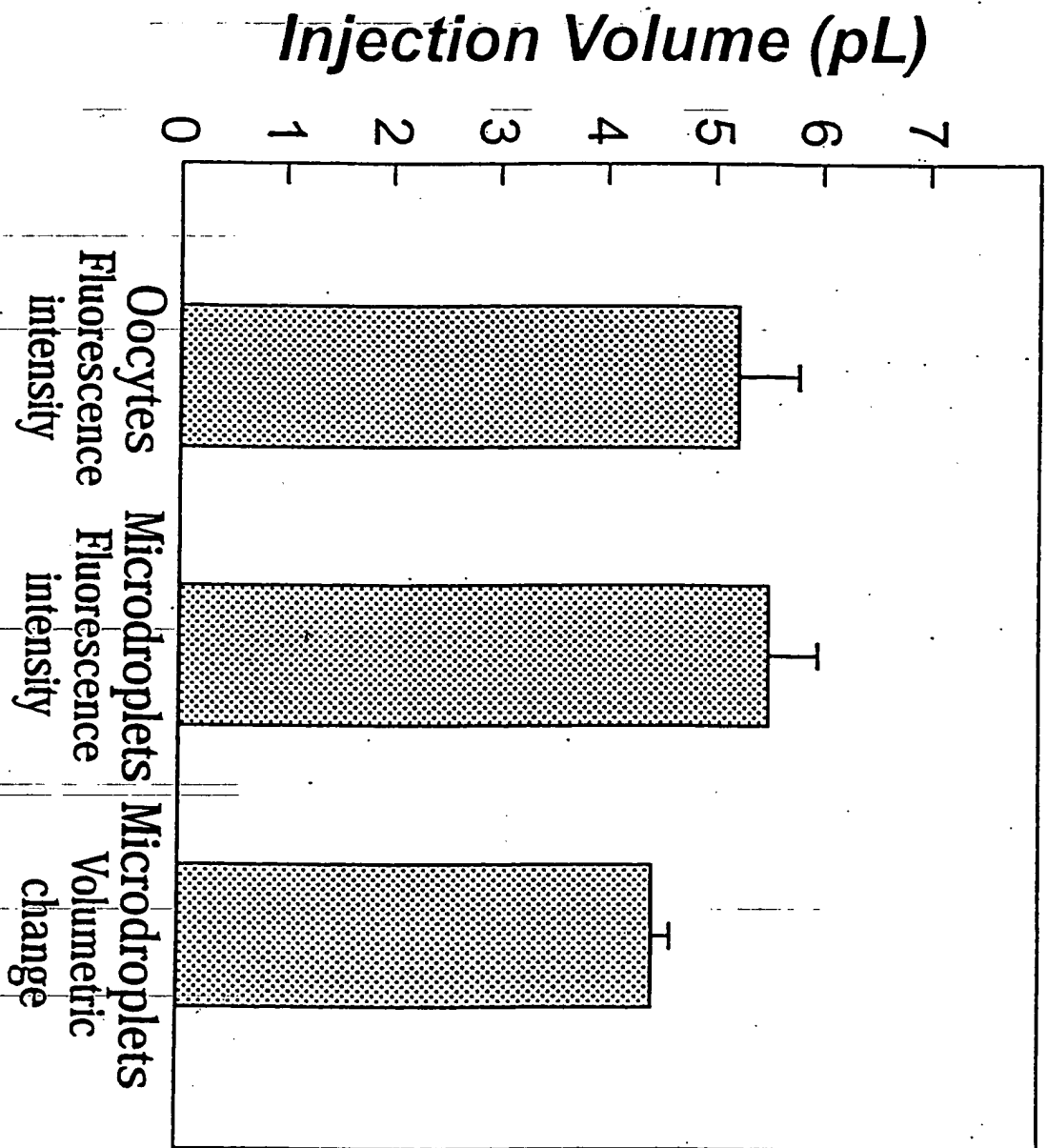


Figure 5

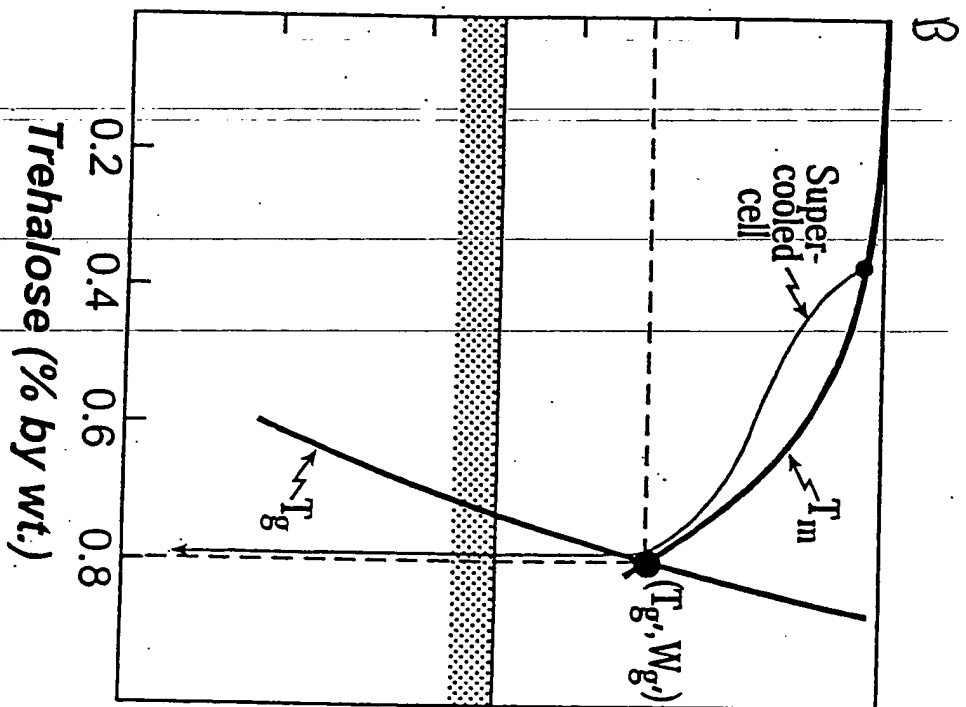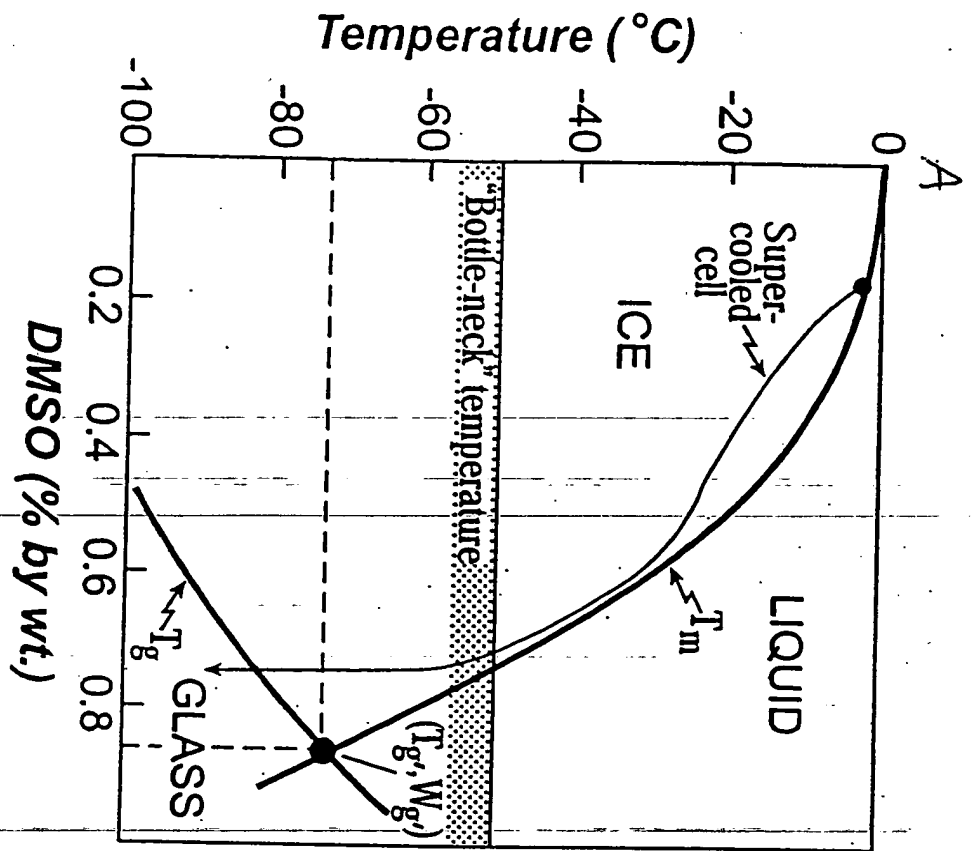


Figures 6A – 6B

Component	Modified HTF		
	HTF	Isotonic	Hypertonic
NaCl (mM)	101.60	↘ 95.00	106.70
KCl (mM)	4.69	↗ 4.78	5.37
KH ₂ PO ₄ (mM)	0.37	↗ 0.38	0.43
MgSO ₄ • 7H ₂ O (mM)	0.20	0.20	0.22
NaHCO ₃ (mM)	25.00	25.00	28.07
CaCl ₂ • 2H ₂ O (mM)	2.04	↘ 2.00	2.25
Lactate (mM)	21.40	↘ 20.00	22.46
Pyruvate (mM)	0.33	0.33	0.37
D-Glucose (mM)	2.78	2.78	3.12
L-Glutamine (mM)	-----	↗ 1.00	1.12
BSA (mg/ml)	4	4.00	4.00
EDTA	-----	↗ 0.01	0.01
Essential amino acids	-----	↗ 0.5x	0.5x
Non-essential amino acids	-----	↗ 0.5x	0.5x
Phenol red (% w/v)	0.001	0.001	0.001
Gentamicin (mg/ml)	50.00	50.00	50.00
Osmolarity (mosm)	285.00	285.00	320.00

Figure 7

Culture Condition	No. of		Percent Two-cell	Percent Blastocyst
	No. of Experiments	Inseminated Oocytes		
<i>Modified HTF, isotonic</i>	5	58	90% (52/58)	87% (51/58)
<i>Modified HTF, hypertonic</i>	9	234	90% (211/234)	86% (202/234)
<i>[Tre]ⁱⁿ = 0.07M</i>	3	32	78% (25/32)	66% (21/32)
<i>[Tre]ⁱⁿ = 0.15M</i>	5	69	67% (46/69)	29% (20/69)

Figure 8

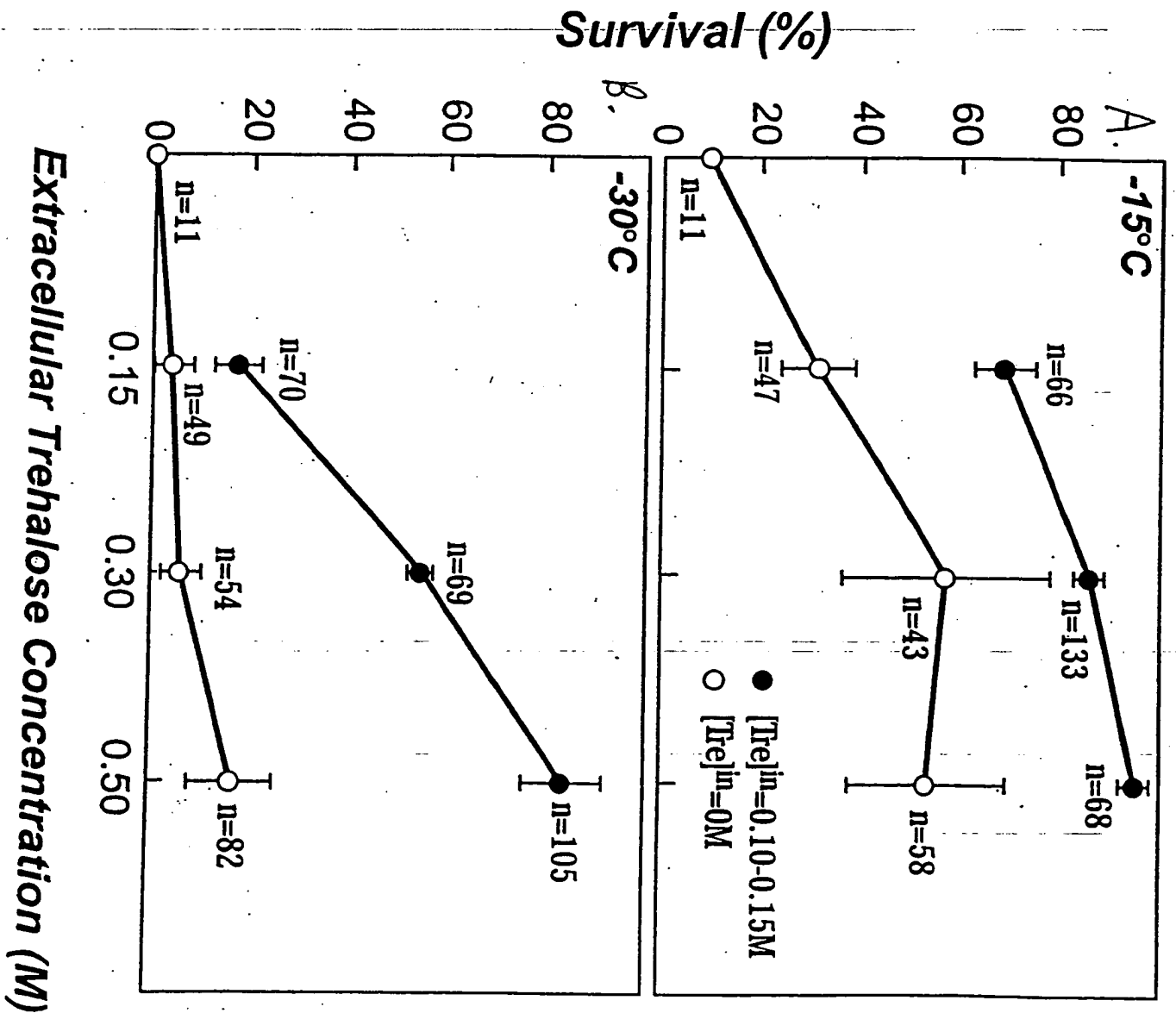

Figures 9A – 9B

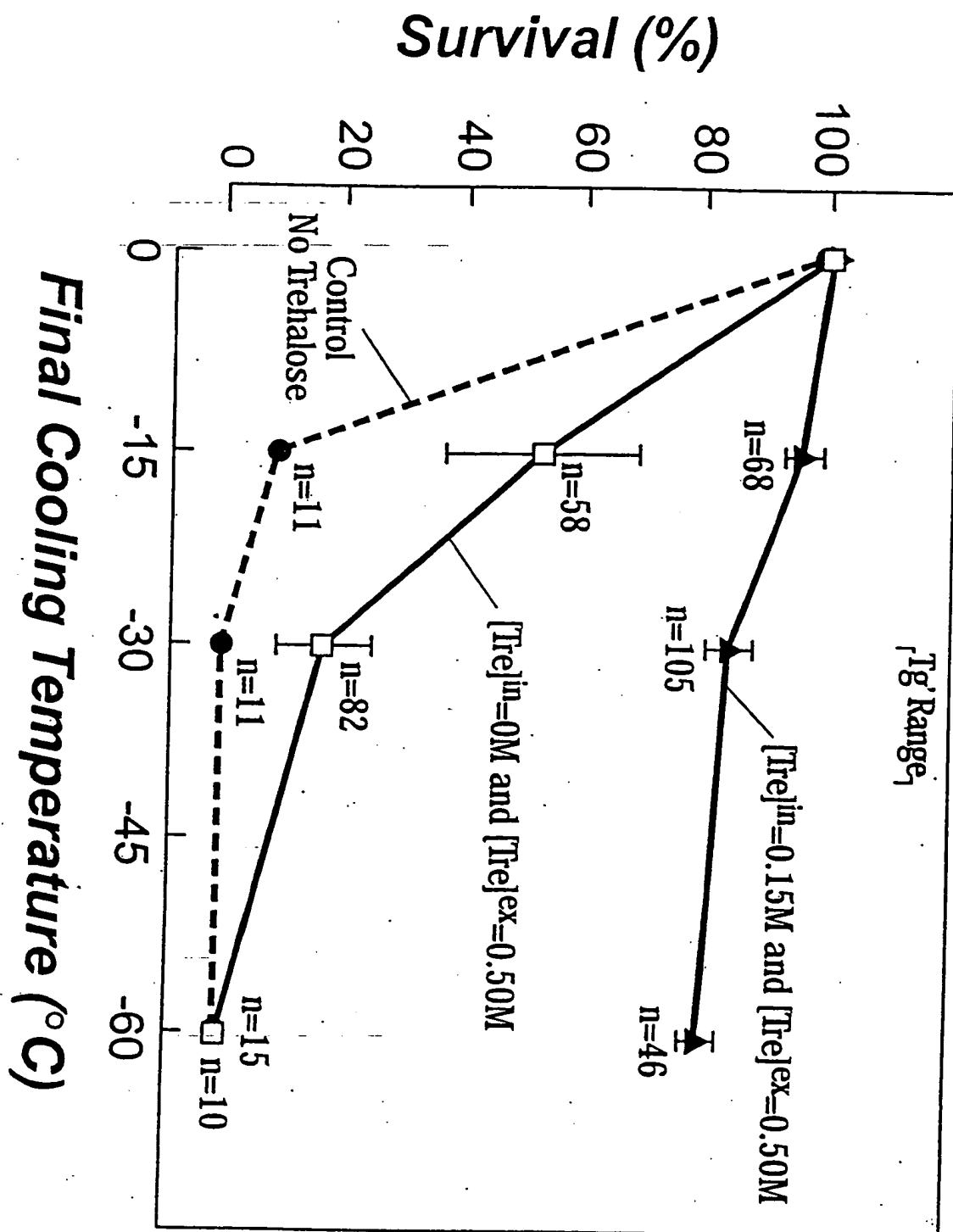


Figure 10



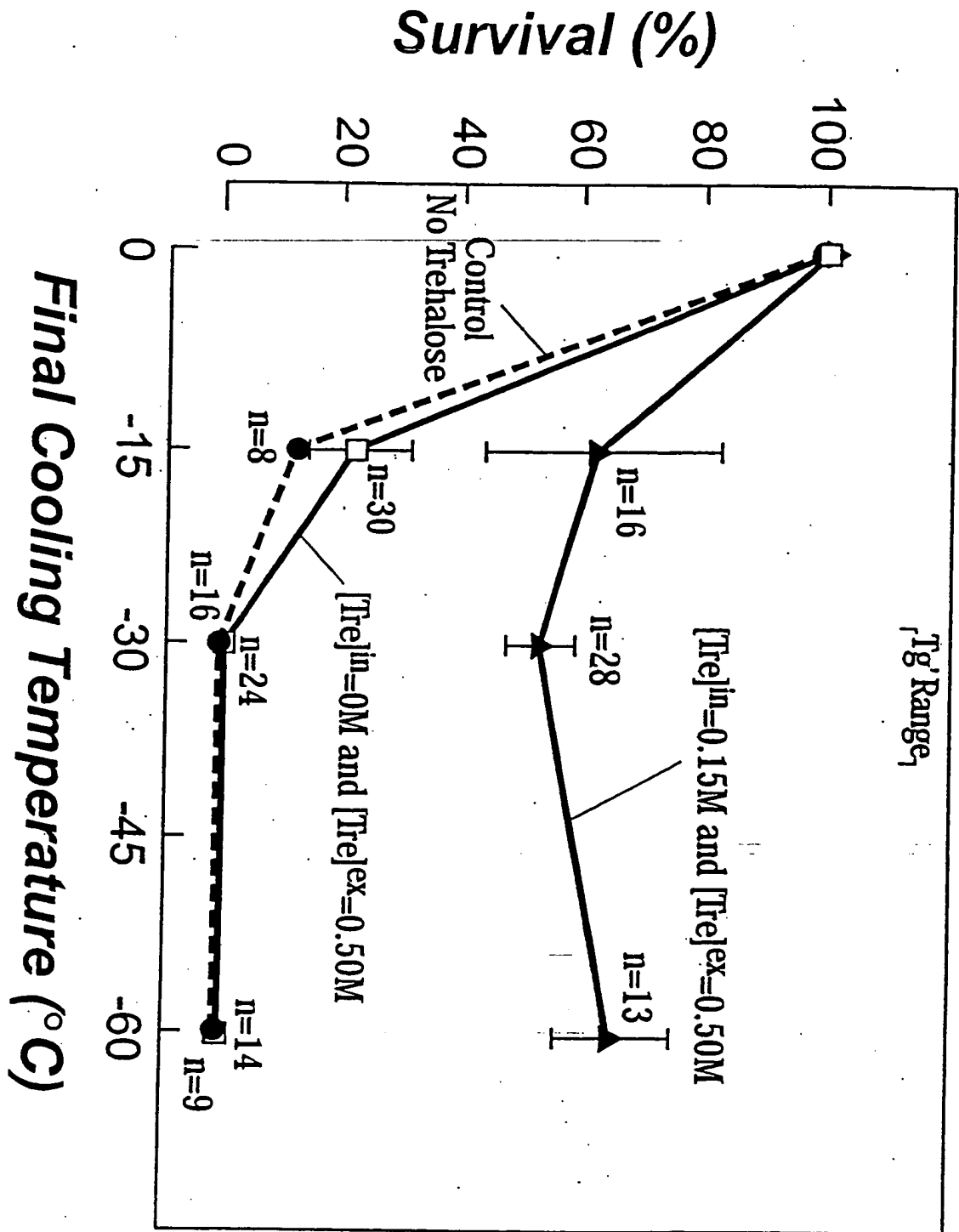


Figure 11

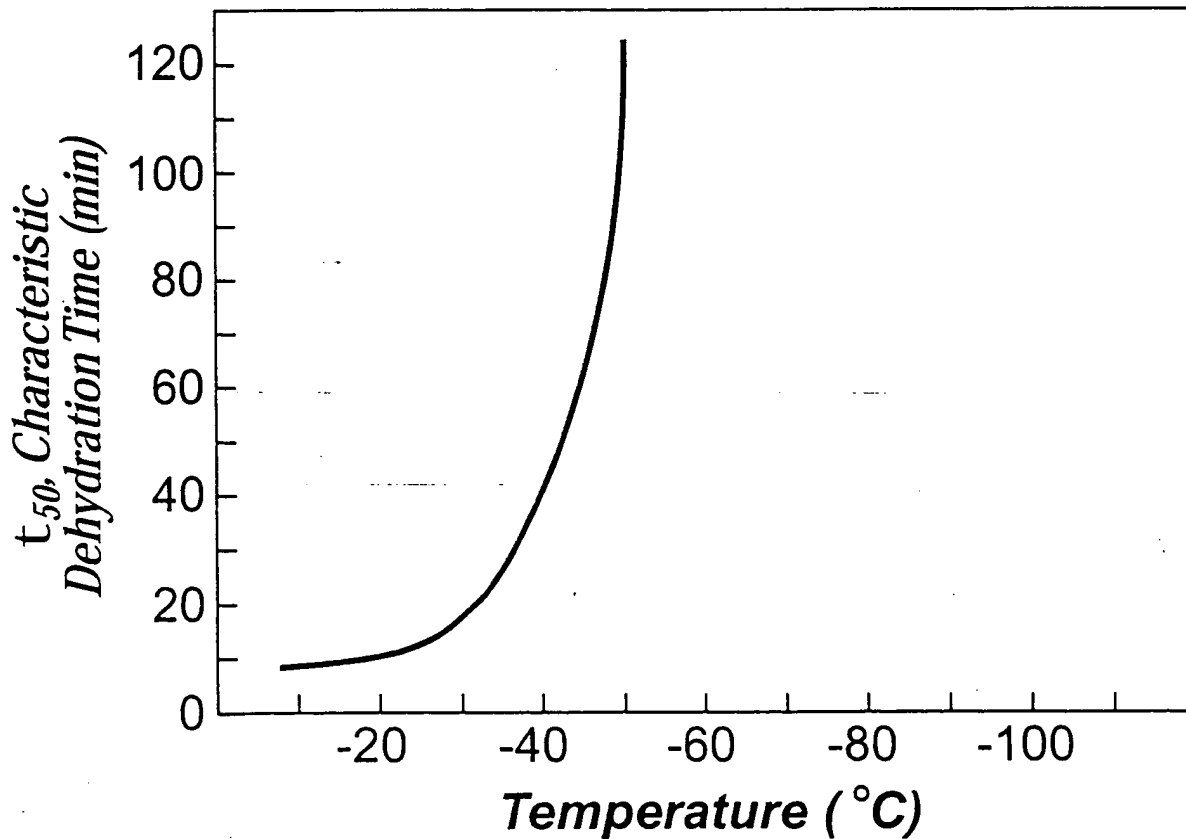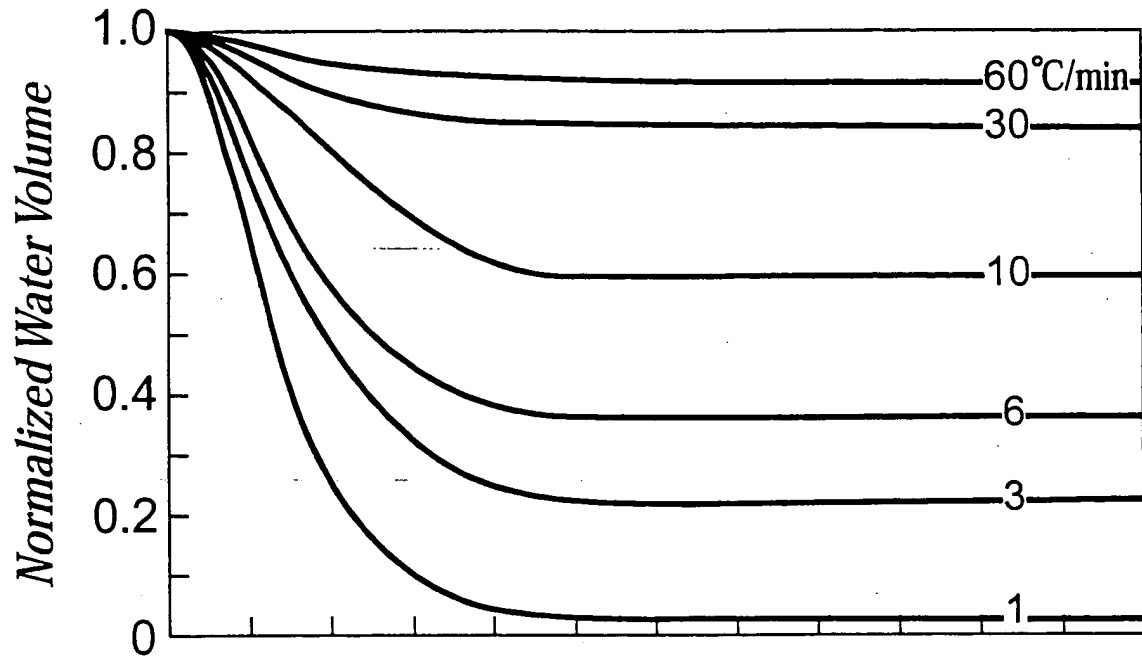

Figures 12A – 12B

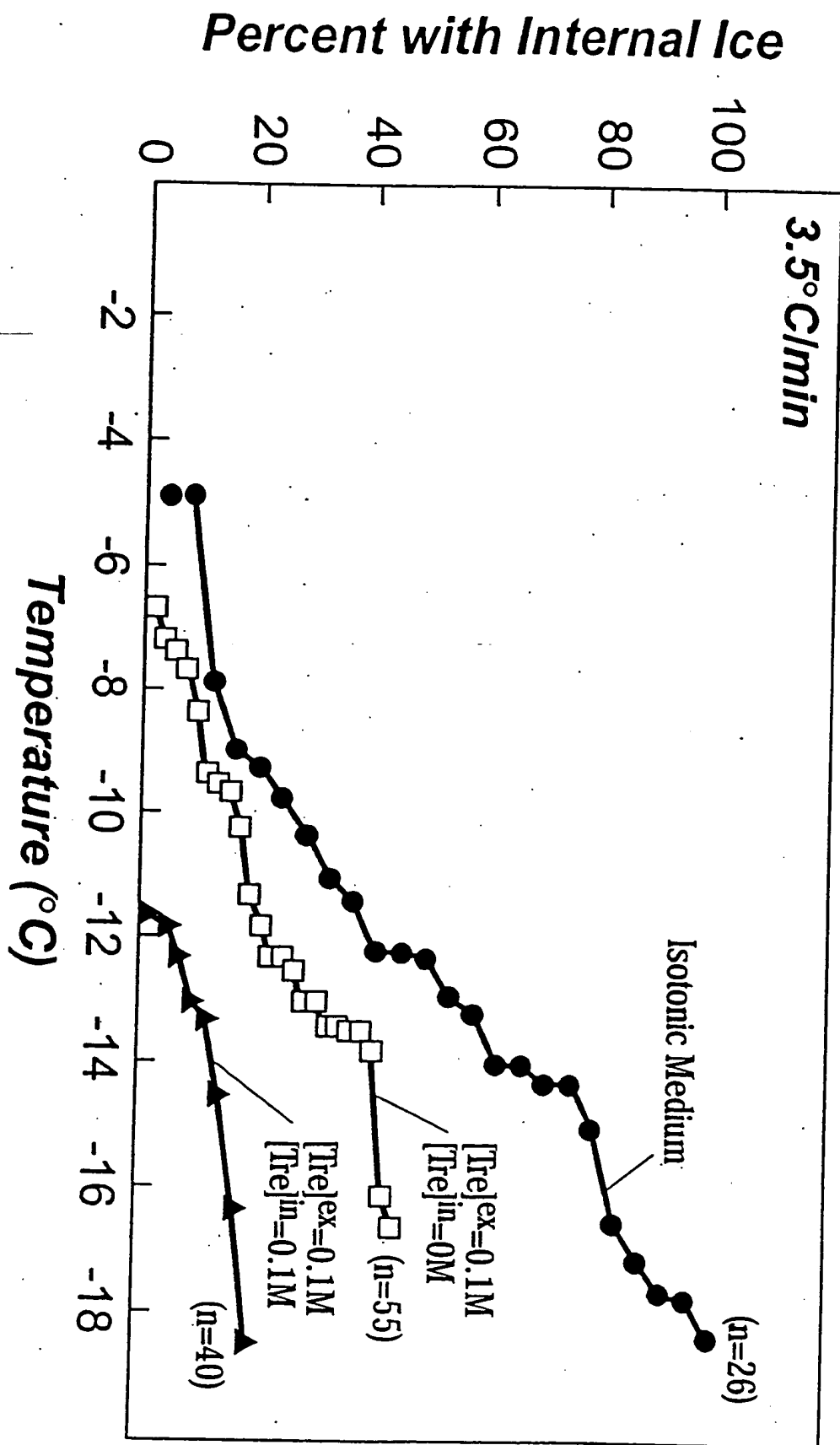


Figure 13

Figure 14

Solute	M _w	T _v °C
erythrose	120.1	-50
threose	120.1	-45.5
erythritol	122.1	-53.5
thyminosse (deoxyribose)	134.1	-52
ribulose	150.1	-50
xylose	150.1	-48
arabinose	150.1	-47.5
lyxose	150.1	-47.5
ribose	150.1	-47
arabitol	152.1	-47
ribitol	152.1	-47
xylitol	152.1	-46.5
methyl riboside	164.2	-53
methyl xyloside	164.2	-49
quinovose (deoxyglucose)	164.2	-43.5
fucose (deoxygalactase)	164.2	-43
rhamnose (deoxymannose)	164.2	-43
talose	180.2	-44
idose	180.2	-44
psicose	180.2	-44
altrose	180.2	-43.5
glucose	180.2	-43
gulose	180.2	-42.5
fructose	180.2	-42
galactose	180.2	-41.5

Title: Microinjection of Cryoprotectants For Preservation of Cells

Applicant(s): Mehmet Toner et al.

Filing Date: December 31, 2003

Page 14 of 24

Customer No.: 21559

Figure 14 Continued

Solute	M _w	T _g °C
allose	180.2	-41.5
sorbose	180.2	-41
mannose	180.2	-41
tagatose	180.2	-40.5
inositol	180.2	-35.5
mannitol	182.2	-40
galactitol	182.2	-39
sorbitol	182.2	-43.5
2-O-methyl fructoside	194.2	-51.5
β-1-O-methyl glucoside	194.2	-47
3-O-methyl glucoside	194.2	-45.5
6-O-methyl galactoside	194.2	-45.5
α-1-O-methyl glucoside	194.2	-44.5
1-O-methyl galactoside	194.2	-44.5
1-O-methyl mannoside	194.2	-43.5
1-O-ethyl glucoside	208.2	-46.5
2-O-ethyl fructoside	208.2	-46.5
1-O-ethyl galactoside	208.2	-45
1-O-ethyl mannoside	208.2	-43.5
glucoheptose	210.2	-37.5
mannoheptulose	210.2	-36.5
glucoheptulose	210.2	-36.5
perseitol (mannoheptitol)	212.2	-32.5
1-O-propyl glucoside	222.2	-43
1-O-propyl galactoside	222.2	-42

Title: Microinjection of Cryoprotectants For Preservation of Cells

Applicant(s): Mehmet Toner et al.

Filing Date: December 31, 2003

Page 15 of 24

Customer No.: 21559

Figure 14 Continued

Solute	M _w	T _g , °C
1- <i>O</i> -propyl mannoside	222.2	-40.5
2,3,4,6- <i>O</i> -methyl glucoside	236.2	-45.5
isomaltulose (palatinose)	342.3	-35.5
nigerose	342.3	-35.5
cellobiulose	342.3	-32.5
isomaltose	342.3	-32.5
sucrose	342.3	-32
gentiobiose	342.3	-31.5
laminaribiose	342.3	-31.5
turanose	342.3	-31
mannobiose	342.3	-30.5
melibiose	342.3	-30.5
lactulose	342.3	-30
maltose	342.3	-29.5
maltulose	342.3	-29.5
trehalose	342.3	-29.5
cellobiose	342.3	-29
lactose	342.3	-28
maltitol	344.3	-34.5
isomaltotriose	504.5	-30.5
panose	504.5	-28
raffinose	504.5	-26.5
maltotriose	504.5	-23.5
nystose	666.6	-26.5
stachyose	666.6	-23.5

Title: Microinjection of Cryoprotectants For Preservation of Cells

Applicant(s): Mehmet Toner et al.

Filing Date: December 31, 2003

Page 16 of 24

Customer No.: 21559

Figure 14 Continued

Solute	M_w	T_g °C
maltotetraose	666.6	-19.5
maltopentaose	828.9	-16.5
α -cyclodextrin	972.9	- 9
maltohexaose	990.9	-14.5
maltoheptaose	1153.0	-13.5

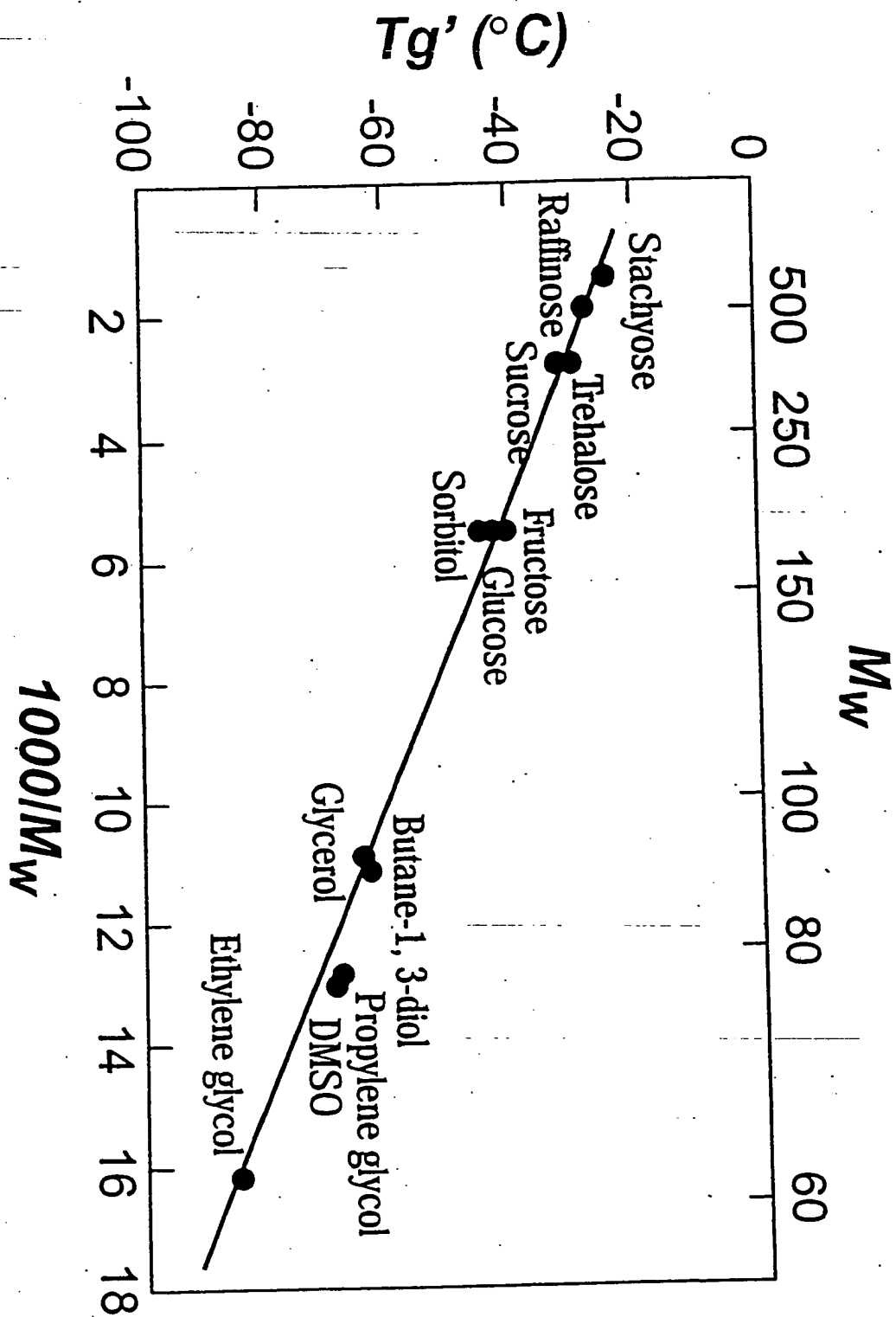


Figure 15

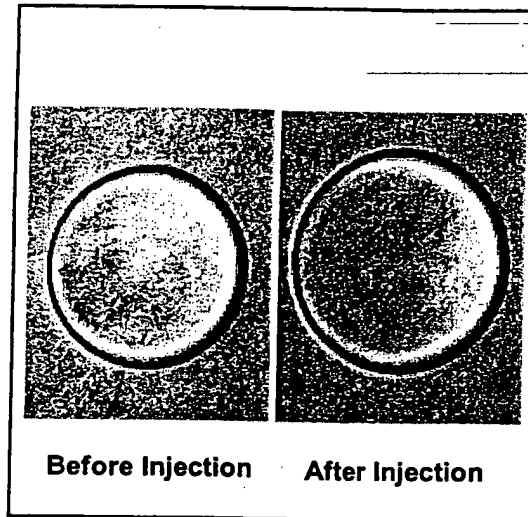


Figure 16

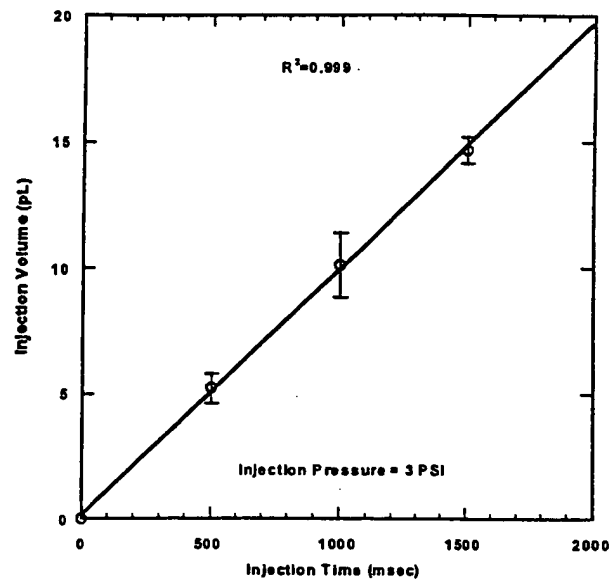


Figure 17

Title: Microinjection of Cryoprotectants For Preservation of Cells

Applicant(s): Mehmet Toner et al.

Filing Date: December 3rd, 2003

Page 21 of 24

Customer No.: 21559

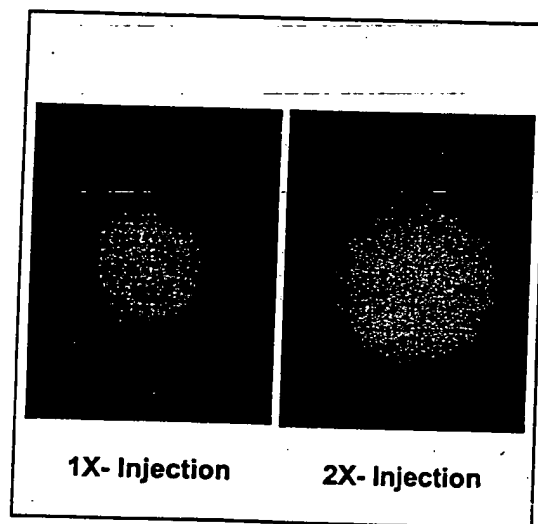


Figure 18

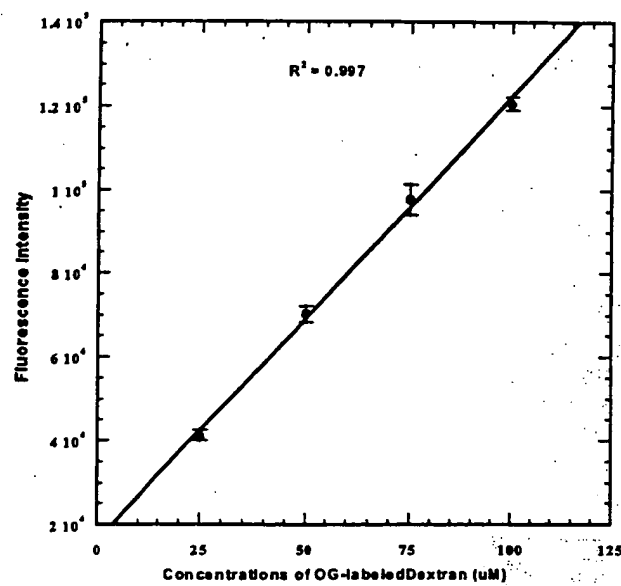


Figure 19

Component	Hyper 320	Hyper 340	Hyper 360
NaCl (mM)	107.02	115.90	123.50
KCl (mM)	5.38	5.83	6.21
KH ₂ PO ₄ (mM)	0.43	0.46	0.49
MgSO ₄ • 7H ₂ O (mM)	0.23	0.24	0.26
NaHCO ₃ (mM)	28.16	30.50	32.50
CaCl ₂ • 2H ₂ O (mM)	2.25	2.44	2.60
Lactate (mM)	22.53	24.40	26.00
Pyruvate (mM)	0.37	0.40	0.43
D- Glucose (mM)	3.13	3.39	3.61
L-Glutamine (mM)	1.13	1.22	1.30
EDTA (mM)	0.01	0.01	0.01
Phenol red (mM)	0.03	0.03	0.03
BSA (mg/mL)	4.00	4.00	4.00
Gentamicin (mg/ml)	50.00	50.00	50.00
Essential amino acids ⁱ (mL/L)	10.00	10.00	10.00
Non-essential amino acids ⁱⁱ (mL/L)	5.00	5.00	5.00
Osmolality	320	340	360

Figure 20

ⁱ [50X Solution, Gibco]

ⁱⁱ [100X Solution, Gibco]

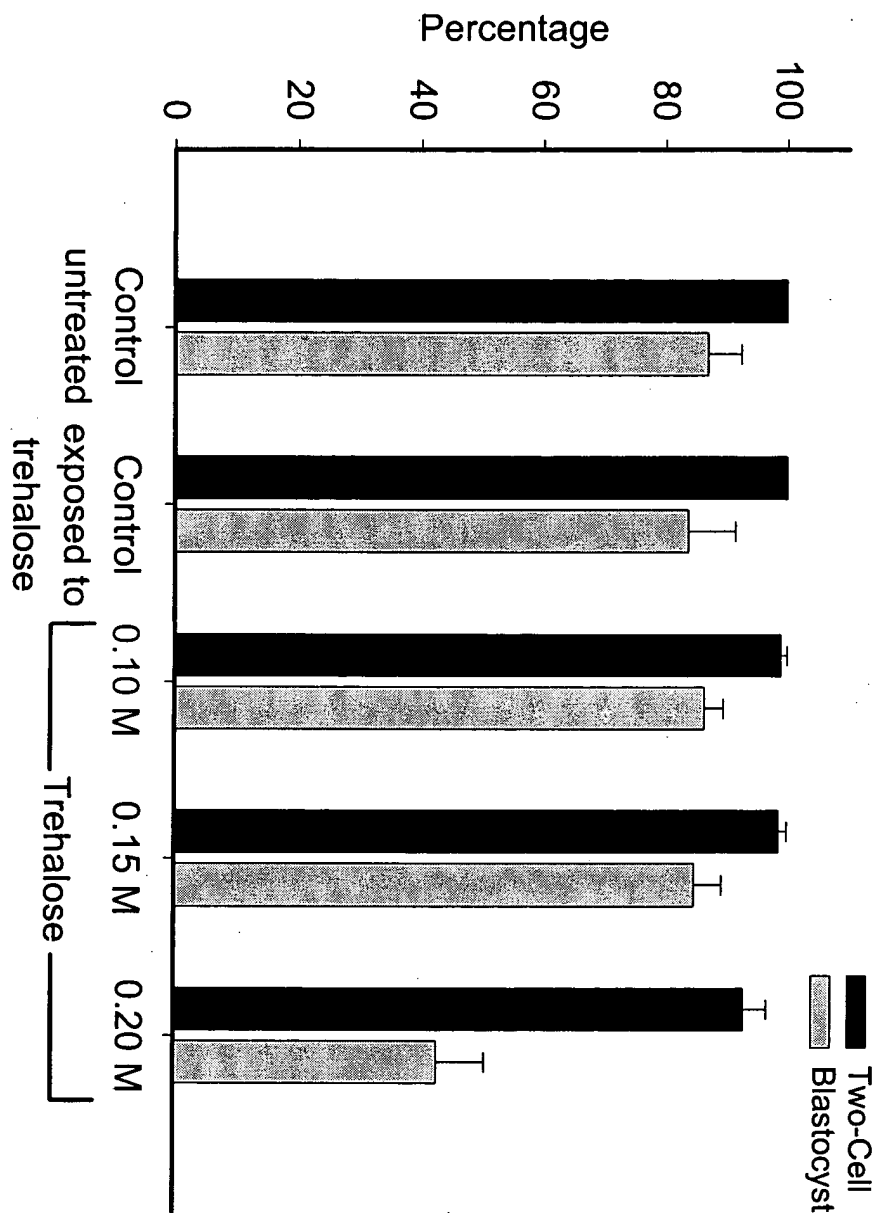


Figure 21